

**"It healed my knees within less than a week. Once he gave me the GlideWear, it was like night and day."**

*Pauline Wiirre*  
Bilateral Transtibial Amputee

**"When I'm using [GlideWear], it totally eliminates my need to take off the limb for periods of rest."**

*LCpl Mar Seifert, USMC, Retired*  
Transtibial Amputee

**"I love the stuff. I can't talk about it enough."**

*Dohn Fadden*  
Transtibial Amputee

**"Simple solutions, like the use of a small GlideWear patch, can allow amputees to gain the full benefits of their advanced prosthetic componentry by allowing them to spend additional time in their prosthesis without the incidence of skin breakdown."**

*Charles W. Kuffel*  
MSM, CPO, FAAOP  
President & Clinical Director,  
Arise Orthotics & Prosthetics



Earn Distance Learning  
Continuing Education Credits


**1.5 Credits**

### CLINICAL EDUCATION CONTENT

Advanced information for healthcare providers and end-users seeking more information about the clinical benefits and specific indications for recommending and using the Tamarack® Prosthetic Liner Patch with GlideWear™ technology.

For online coursework and accredited post-presentation quiz, please visit:

[beckerorthopedic.com/Education/ContinuingEducation](http://beckerorthopedic.com/Education/ContinuingEducation)

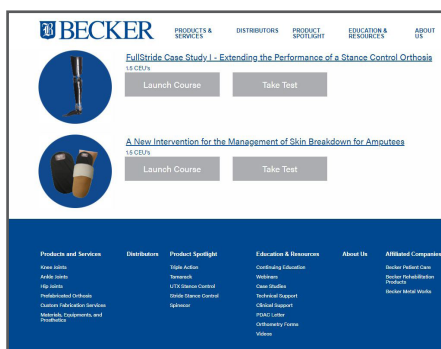

**GlideWear™**  
TECHNOLOGY

## Tamarack® Prosthetic Liner Patch

for Immediate Relief and Lasting Protection  
from Skin Breakdown



**Tamarack®**  
HABILITATION TECHNOLOGIES, INC.

MADE IN MINNESOTA ★ U.S.A. ★

Tamarack Habilitation Technologies, Inc.  
1670-94th Lane NE, Blaine, MN 55449  
Tel: 866.795.0057

E-mail: [info@tamarackhti.com](mailto:info@tamarackhti.com)

ISO13485: 2016 Registered QMS

© 2023 Tamarack Habilitation Technologies, Inc.  
GlideWear is a registered trademark of MIPS, manufactured in the USA by Tamarack HTI.  
Patented US8646459 - US9241853 - GB2478206 - CA2738162 - Other Patents Pending. All Rights Reserved.

## FEATURES & BENEFITS

### Easy to Use

Place the Tamarack® Prosthetic Liner Patch directly against the area of concern, before donning liner over the top of the patch and limb. Use a clean patch daily to protect against future breakdown.

### Stays in Place All Day Long

Tamarack® Prosthetic Liner Patches are held in place by the gel liner or prosthetic socket; ensuring uninterrupted skin protection all day long.

### Reduces Shear

Patented dual-layer GlideWear™ technology absorbs friction and related shear forces at the skin - prosthetic liner interface.

### Durable & Reusable

Machine or hand wash frequently to keep clean and maintain effectiveness.

### Everyone Saves

Amputees save money & energy otherwise spent purchasing and applying messy creams, powders and bandages.

Healthcare professionals reduce non-billable time and office expense associated with socket modifications and wound care efforts.

## SPECIFICATIONS

**Content:** 85% Nylon, 15% Spandex.  
Does not contain latex.

### CARE INSTRUCTIONS

Wash with mild detergent.



### Available Sizes

#### GW-PRO-2SB

##### Small Oval

2.5 in x 4.75 in  
(6 cm x 12 cm)

#### GW-PRO-2LB

##### Large Oval

4.5 in x 7.5 in  
(11 cm x 19 cm)

#### GW-PRO-2CB

##### Circle

4 in diameter  
(10 cm diameter)



## INDICATIONS

The Tamarack® Prosthetic Liner Patch with GlideWear™ technology is **indicated for lower limb amputees who experience painful skin irritation and breakdown** on their residual limb from daily use of a prosthesis.

When used in conjunction with a prosthetic liner, GlideWear™ provides **immediate relief and lasting protection** from skin breakdown caused by friction & shear forces between skin tissue, bony prominences, and the prosthesis.

## APPLICATION & PLACEMENT TIPS

- Contact your Prosthetist to determine if the Tamarack® Prosthetic Liner Patch with GlideWear™ technology is appropriate for your skin care needs.
- **Place the Tamarack® Prosthetic Liner Patch between the prosthetic liner and skin** - adjacent to the area of concern - while donning the liner.
- **Do not place directly over an open wound.**
- **For placement in distal applications**, lay the Tamarack® Prosthetic Liner Patch over the rolled edge of the liner before donning the liner onto the limb.
- **Use a clean and dry patch daily** to ensure lasting protection from discomfort, redness, and skin breakdown.



BEFORE

FIBULAR HEAD / HAMSTRING APPLICATION

AFTER



For more information, visit

[www.tamarackhti.com](http://www.tamarackhti.com)

BEFORE

AFTER

ANTERIOR DISTAL TIBIA APPLICATION



BEFORE

AFTER

PATELLA APPLICATION

