

VACUUM INTEGRATED PUMP

VIP



VIP TUTORIALS



DESCRIPTION

The ALPS™ VIP is an air expulsion integrated pump for elevated vacuum for prosthetic socket. The use of VIP provides patients continuous total contact when donning the prosthesis. This is essential for an even weight distribution and to distribute the suspension of the prosthesis over the whole surface of the residual limb. Elevated vacuum created by the VIP stabilizes the residual limb, overcoming all issues connected with the residual limb changing volume during the day.

ADDED BENEFITS



Efficiently deals with stump volume fluctuation which in turn leads to improved prosthesis control and comfort



Reduces the intensity of peaks in pressure on the patient's limbs, improving comfort and providing greater freedom of movement



Helps preserve residual limb skin health by improving circulation and consequently increasing limb oxygenation



Does not affect the cosmetic appearance of the prosthesis

Easy to
install

Improves
blood flow

Improves
comfort

Improves
walking

Quiet

Light



#makinglivesbetter

MAKING LIVES BETTER

AVAILABLE IN 2 DIFFERENT VACUUM LEVELS

HIGH VACUUM (Red)
10-13 in of mercury



We highly recommend using the **mini valve** to prevent distal distraction

LOW VACUUM (Green)
5-8 in of mercury



**CHOOSE LEVEL
BASED ON:**




ACTIVITY LEVEL




SKIN CONDITION

WITHOUT MOUNTING PLATE

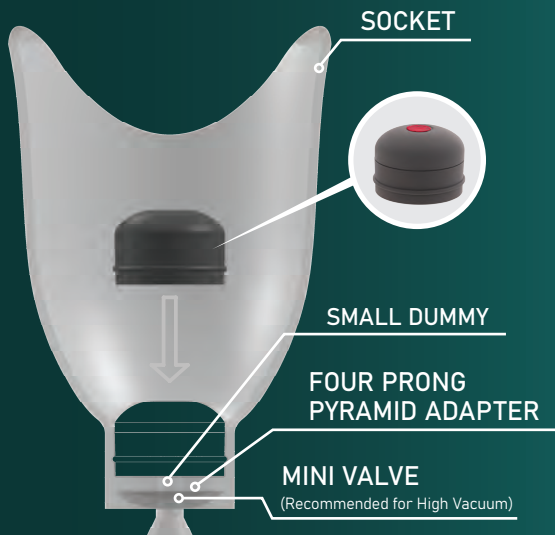

WITH MOUNTING PLATE

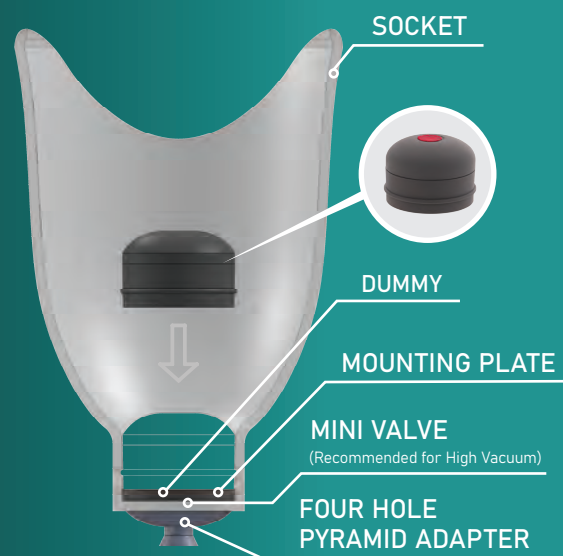

PRODUCT SPECS

	HIGH VACUUM (RED)	LOW VACUUM (GREEN)
PART #	VIP - DH	VIP - DL
DIAMETER	66MM	66MM
HEIGHT	48MM	48MM
WEIGHT	75G	64G
MATERIAL	SILICONE	SILICONE

USE CODES BELOW WHEN ORDERING

PART NAME	PART #
DUMMY	VE319 - D
SMALL DUMMY	VE319 - SD
MOUNTING PLATE	VE819 - P
LOCK (SMALL)	VIP-L1120-S
MINI VALVE (Highly Recommended)	MV-009
MINI VALVE KIT	LMV-009
PIN (1.9 IN)	VIP-92L
PIN (1.2 IN)	VIP-92S



T: 727.528.8566 | 800.574.5426 | F: 727.528.8862

www.easyliner.com | info@easyliner.com

2895 42nd Avenue North St. Petersburg, Florida 33714

